INSTRUCTIONS FOR USE all Food:

This sampling kit is a send-in kit for the measurement of food specific IgG antibodies in human **capillary sample**. The sample is analysed in vitro in a specialized diagnostic laboratory.

Content of sample sampling kit:

Material for sampling (sample sampling pouch)

- 1 alcohol wipe
- 1 band aid
- 1 swab
- 1 collection tube
- 1 transport tube
- 2 lancets

1 Strip with two Test ID / QR code labels

Written Material

1 Patient form

1 Instructions for use

Please note:

- Please read the Instructions carefully before you begin.
- If you have difficulties drawing sample or if the kit is incomplete, you can contact our customer service team with complete confidence.

Warnings:

- Keep the kit out of reach of children.
- kit may only be used by adults over the age of 18.
- If you have a congenital or acquired sample clotting disorder or take anticoagulant medication, use the kit only after consulting your doctor.
- The kit and especially the sterile lancets are exclusively for your own use and may only be used once.
- If you have problems regarding the sight of or the taking of sample, take the sample only in the presence of another person.
- Any serious incident that has occurred in relation to the use of this kit shall be reported to the manufacturer and the competent authority of the State in which you are living.

Do not use the kit, if you have a sample-borne disease, e.g. Hepatitis or HIV!

Instructions:


1. [to be done: picture of patient form]

Before you start, you need to fill in the patient form. Your test kit contains two stickers with your unique test-ID and a QR code. Place one sticker on the patient form.

2. If you have cold hands, wash your hands in hot water for about a minute. Dry your hands. Additionally, you may swing your arm in circles for 10 times. Both actions will increase sample circulation in your fingers.


3. Wipe your middle or ring finger to be punctured with the alcohol wipe. Wait a few seconds until the finger is completely dry.

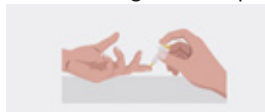
4. Remove the protective cover of the lancet by rotating and pulling.

CAUTION: Do not use a lancet with a loose or missing protective cover.

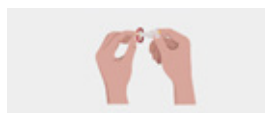5. Place your hand on a flat surface. Put the lancet to the side of your finger tip and press the lancet firmly against the finger, until the lancet needle is released. **CAUTION:** the lancet will only release once.

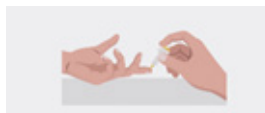
6. We recommend that you stand up at this stage to let the sample in your arm flow downward. Help the sample drops come out by massaging the finger toward the needle puncture.

7. Open the collection tube by removing the cap. Transfer the sample drops into the test tube until the sample reaches the marked line. Standing up will help with sample flow. Close the cap after finishing.


8. It is **very important** that you fill up the tube to the marked line, otherwise the lab will not be able to complete the test. Use the extra lancet if you need to puncture another finger.

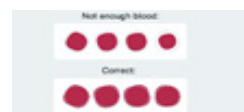9. Place the sticker with your unique test-ID on the collection tube. Place the collection tube inside the transport tube and close it.

10. When the sample has dried for 2 hours, put the paper in the plastic cover it arrived in.

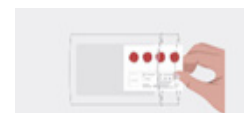
11. Using a return envelope supplied by your sales agency, the sample will be send to the laboratory.

Icon explanation:

Consult instructions for use.

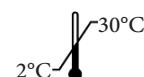Storage between 2°C - 30°C.

Use by.

Caution.

Article number.

Do not reuse

Batch.

Manufacturer.

Content sufficient for <1> test.

In-vitro diagnostics.

CE Marking of conformity.

GA Generic Assays GmbH,
Ludwig-Erhard-Ring 3, 15827
Blankenfelde-Mahlow OT
Dahlewitz, Germany.

GA Generic Assays GmbH
Ludwig-Erhard-Ring 3
15827 Blankenfelde Mahlow / OT Dahlewitz
GERMANY
t: +49 33708 9286 0 - +49 1723993687
info@genericassays.com
www.genericassays.com
Trade registry HRB 16038
District court Potsdam

Medipan GmbH
Ludwig-Erhard-Ring 3
15827 Blankenfelde Mahlow / OT Dahlewitz
GERMANY
t: +49 33708 4417 32 - t: +49 33708 4417 25
s.landau@medipan.de
www.medipan.de
Trade registry HRB 6195
District court Potsdam

King Abdulla District Building no. 8178 East
ring road sub. 12424 -4352
Email: info@hmc.health
Saudi Arabia
Phone: + 966112416899
Fax: +966112708255
www.hmc.health

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